


INTRODUCING THE PARALLEL™


The first disposable nose clip that
FITS and STAYS on EVERY nose.

Quality is so good, you think it's a permanent nose clip.

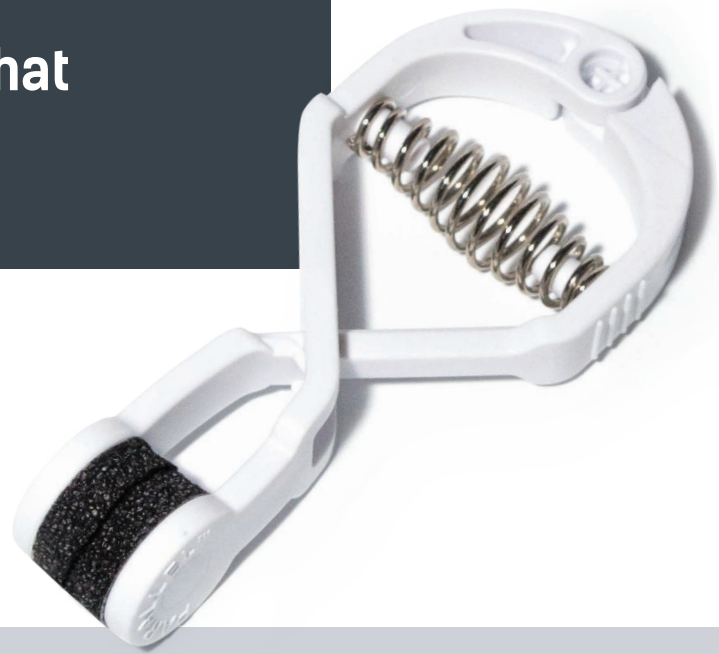

THE PARALLEL™ fixes the two major complaints about nose clips:

1. Fitting every patient's nose
2. Staying on every nose

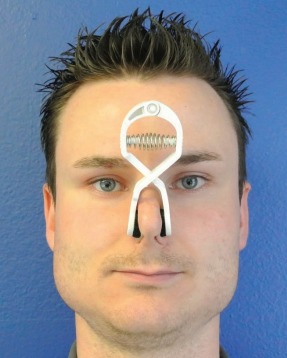

How is the PARALLEL different:

FITS ALL NOSE SIZES, EVEN EXTRA-WIDE

The standard nose clip can only spread to about 3 cm, much smaller than many noses. The Parallel spreads to over 5.5 cm (an inch more) and can capture the noses you now hold with your hand.

THE PARALLEL is the perfect nose clip for small children through the largest adult.

STAYS ON, SLIP-PROOF

Standard disposable nose clips angle wide as you spread them. When placed on most noses, the excessive angle means the clip slips right off your patient's nose.

THE PARALLEL opens in nearly a parallel motion which minimizes the angle of the spread. Less angle, no slipping.



Standard nose clip's narrow opening doesn't fit all noses.

The PARALLEL's wider opening fits ALL nose sizes!



Standard nose clip is angled and slips off the nose.

The PARALLEL's design stays on the nose.

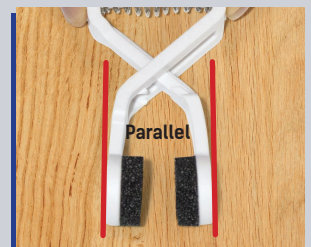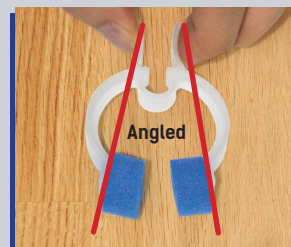

OTHER FEATURES

- Designed to block air leaking from the nose
- THE PARALLEL's metal spring seals the nostrils with consistent force
- Provides unobstructed visibility to screens and equipment
- The textured foam pads prevent slipping on the skin, while being extremely comfortable
- Pads meet ISO 10993-1 biocompatibility requirements



MANUFACTURED FOR:
12th Man Technologies, Inc.
7245 Garden Grove Blvd., Suite C
Garden Grove, CA 92841
1.714.705.4585 • www.12thmantec.com

PN 50-50043 Rev. A • Made in China • © 2016 12th Man Technologies, Inc.

PRODUCT INFORMATION:

Disposable Nose Clips – Box of 100 (PN 50-90059)

ALSO AVAILABLE:

Parallel-R™ with replaceable Nose Pads

- Box of 100 (10 pairs Nose Pads with each Clip) (PN 50-90074)
- 100 Replacement Nose Pads (50 pairs) (PN 50-90076)

For More Information and to Order:



Novus Medical Inc.
1-866-926-9977
www.novusmedical.ca
info@novusmedical.ca

